


## Selected Electrode Placement Chart

The sticky electrode pads need to be placed in the correct position for maximum pain relief. It is preferable to obtain assistance from either your birthing partner or medical adviser when positioning the pads. Please read the following instructions carefully:

- Carefully peel all of the pads from the plastic backing using your forefinger and thumb. Do not pull on the lead wires.
- Position the top pads (connected to the mauve coloured leads) approximately 4cm either side of the spine and 4-5cm above the imaginary navel line as seen in the diagram.

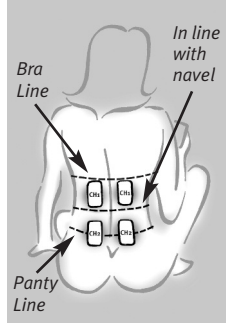

Reference: Cluett E. (1994) SRN RM ADM PGCEA. Analgesia in Labour: A review of the TENS method. Professional Care of Mother and Child. Mar/4 (2): 50-2.

7


## Useful Information

- Continue stimulation for as long as necessary, unless you experience discomfort.
- If you stimulate for a number of hours ensure the pads do not become dry. If this occurs, switch off the unit, disconnect the electrodes, and apply a small drop of water to them.
- Always ensure that the unit is switched off before you remove all the leads and pads. Do not pull the wires. ALWAYS remove the jack plugs from their sockets by holding the plugs between the forefinger and thumb.

10

- Position the bottom two pads (connected to the white coloured leads) approximately 4cm either side of the spine and 4-5cm below the imaginary navel line as seen in the diagram. NB. It makes no difference whatsoever whether the mauve pads are at the top and the white ones are below or vice versa, as long as they are in horizontal pairs.

- The Obi TENS can also be used as a standard TENS unit for relief of numerous conditions, including after pains and post operative pains. It is worth bearing this in mind if you have a caesarean section. Please seek medical advice first.

 Before first use of the Obi TENS, please familiarise yourself with the precautions and contraindications contained on page 1 and 2.

## Using Your Obi TENS





- You are now ready to start using your Obi TENS. Switch the unit on by pressing the on/off button. The Obi TENS will always start up in BURST mode.
- Begin by pressing the Intensity Up key. Each time you press this key, a small 'half-moon' shape will appear on the screen increasing in number and size as the intensity is increased. You will begin to feel a pulsating sensation which will burst in time with the 'half-moon' shape as shown on the LCD. You will feel this sensation in both sets of pads.
- BURST (use between contractions):** This is the first mode that you will use. It feels like a bursting,

8

- tingling sensation i.e. the pulses are off and on periodically. This mode will help to promote those all important natural pain-killing chemicals known as 'endorphins' and 'encephalins'.
- BOOST (use during contractions):** When experiencing a contraction, simply press the 'BOOST' key. The Obi TENS will then switch automatically into the BOOST mode for that extra surge of power required to combat the pain during contractions. This will feel like a continuous sensation.
- Once the contraction has passed, simply press the 'BOOST' key. The Obi TENS will return to the original BURST mode. Continue this cycle during your early stages of labour.

KEY: ■ ■ = first pair of electrodes □ □ = second pair of electrodes.

N.B. Electrodes supplied with this unit are all the same colour - the black and white squares show suggested electrode placement only.

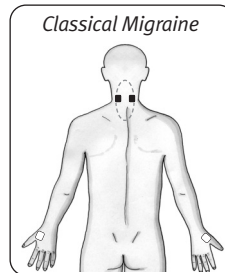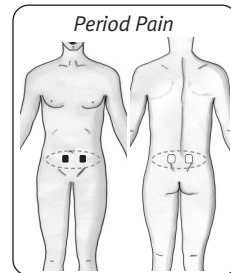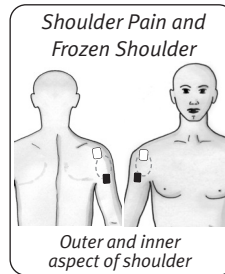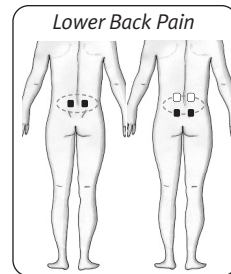

11

## Frequently Asked Questions

- When shall I start using my Obi TENS?**
  - Use as early as you can after the onset of labour to allow time for your body's pain-killing chemicals to rise.
- Can I try the Obi TENS before going into labour?**
  - Yes, place 2 pads on your forearm to feel the different pulse sensations following the instructions on page 6.
- Can I combine the Obi TENS with other medication?**
  - Yes, you can still use analgesics such as gas and air (entinox) or pethidine.
- Can the Obi TENS be used in hospital?**
  - Yes. However, a midwife may ask you to switch it off temporarily.
- Can I use the Obi TENS if I am considering a water birth?**
  - Yes, but not in the water itself.

9

## Troubleshooting Questions

- Why does the stimulation not appear to be as strong after you have used the Obi TENS for a while?**
  - Increase the intensity, you may have become used to a lower setting.
- Why does the stimulation feel strong but ineffective?**
  - You may need to switch off and reposition the electrode pads.
- Why does the stimulation not feel sufficiently strong even on a high setting?**
  - The batteries may need replacing.
- What do I do if there is no stimulation, intermittent stimulation, or weak stimulation felt even with a new battery?**
  - Check electrode contact. Check leads. Check batteries. Phone Babycare TENS on 0845 454 230 9737.



Instructions for use  
**Obi TENS**



## Technical Data

<b>Display</b>	The LCD has a total of 13 display segments. The array of Level Bars consists of 10 LCD Level Bar. In addition, 2 individual LCD segments respectively indicate the Modes (i.e. Burst and Boost). The last LCD segment is used to illuminate the low battery indicator.		
<b>Mode</b>	2 operation modes		
	Pulse (Hz) Frequency	Pulse Width (µs)	Pattern
<b>Burst</b>	32	200	2 bursts per second and 8 pulses per burst
<b>Boost</b>	80	100	Constant pulse
<b>Output Channel</b>	ONE Channel		
<b>Output Voltage</b>	0 – 100V (zero-to-peak) across a 500Ω load; adjustable in 10 steps		
<b>Output Intensity</b>	200mA (max.) into a 500Ω load		
<b>Output Waveform</b>	Symmetrical Mono-phasic rectangular		
<b>Output Jacks</b>	CE Touch-proof jacks		
<b>Battery</b>	2 X1.5V AAA Battery (Type LR03)		
<b>Battery Low Detect</b>	2.2V±0.1V		
<b>Tolerances</b>	All output parameters are subjected to a ±10% tolerance unless other specified		

12



ISO 9001:2008 ISO 13485:2003  
EC DIRECTIVE 93/42/EEC Annex V

**Babycare TENS**

**Babycare TENS**

108 George Lane, South Woodford,  
London E18 1AD

Tel (+44) (0)20 8532 9595

Fax (+44) (0)20 8532 9551

www.babycaretens.com

e-mail customercare@babycaretens.com

Babycare TENS is part of the  
Body Clock group of companies

body clock

CE

0120

August 2012 (5)

recycle

This manual is recyclable

**The Obi TENS is an ultra simple yet powerful obstetric TENS unit perfect for drug-free pain relief during labour. It has 2 pre-set modes; one for use during contractions, the other for use between contractions.**

TENS is recognised as a safe and highly effective method of pain relief and is regularly recommended by medical professionals.

TENS is a drug free method of pain relief, with no known side effects, which can also be used in conjunction with any additional medication if required.

#### Explanation Of Symbols On Unit



Equipment providing a particular degree of protection against electric shock particularly regarding allowable leakage currents having a BF type (floating) applied part.




Warning – refer to accompanying documents i.e. these instructions.



Denotes a product which must be disposed of safely.

### Controls On The Obi TENS

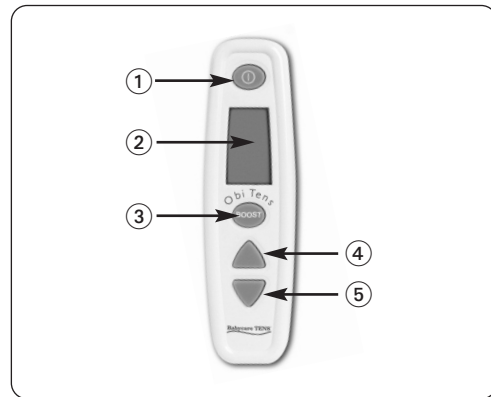


- 1 **POWER** Key – Turns the Obi TENS On and Off.
- 2 **User display** – Indicates the mode in use and the intensity level. (also see diagram below).
- 3 **BOOST** Key –Used for extra surge during contractions.
- 4 **Intensity Up** – Increases the intensity of the pulses.
- 5 **Intensity Down** – Decreases the intensity of the pulses.

### Contents

How Will The Obi TENS Help Your Pain?.....	1
Precautions And Contraindications .....	1
The Obi TENS Controls .....	3
Additional Functions.....	4
Setting Up Your Obi TENS.....	5
Suggested Pad Placement Chart .....	7
Using Your Obi TENS.....	8
Frequently Asked Questions .....	9
Troubleshooting Questions .....	9
Useful Information .....	10
Post-Natal Pad Placement Diagrams .....	11
Technical Data.....	12

### User Display

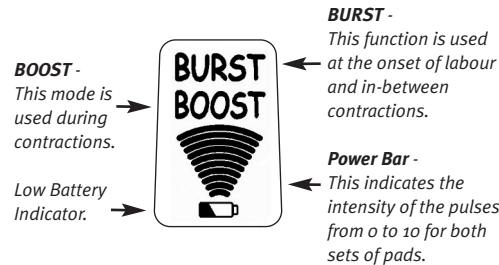

### Additional Functions

- Low Battery Indicator -**  
A flashing battery symbol will be shown on the display.
- Automatic Power Off -**  
When the unit intensity levels are zero, and the Obi TENS has not been in use for about 5 minutes, the unit will shut off automatically.

### How Will The Obi TENS Help Your Pain?

TENS (Transcutaneous Electrical Nerve Stimulation) works by passing mild electrical impulses through the skin, via electrode pads, into the nerve fibres which lie below. The TENS impulses help your body produce its own pain killing chemicals, such as 'endorphins' and 'enkephalins'.

The level of pain relief obtained from the Obi TENS varies from person to person. Some find that the Obi TENS provides all the pain relief required during the whole of labour, whereas others find that other analgesics are also needed during the later stages of labour.



### Precautions And contraindications

1. Do **NOT** use the Obi TENS if you have a cardiac pacemaker.
2. Check with your medical adviser before using if you suffer from a heart condition, epilepsy, undiagnosed pain, have any metal implants, or any doubts whatsoever.
3. TENS must **NOT** be used before the 37th week of pregnancy.
4. Never use TENS to mask undiagnosed pain since this could require urgent treatment.
5. If you are pregnant, do **NOT** place the electrodes over your abdomen.

### Setting Up Your Obi TENS

Check you have the following contents:

- 1 x Obi TENS unit
- 1 x Pack of 4 self adhesive electrodes (40mm x 100mm)
- 1 x Leadwires
- 2 x AAA batteries
- 1 x Carrying pouch
- 1 x Easy-release neck cord
- 1 x Instruction manual

### The Assembly Stage

1. Slide open the battery cover on the back of the unit to expose the battery compartment.
  2. Insert 2 x AAA batteries according to the diagram found in the battery compartment. Replace the cover. Please ensure your unit is switched off.
- a) Take one of the mauve and white leadwires and insert the white jack plug into the socket at the bottom of the Obi TENS unit.

### General Precautions

1. Do **NOT** use this unit without first reading these instructions.
2. Do **NOT** immerse the Obi TENS in any liquid.
3. Do **NOT** place it close to any source of excessive heat or operate it in the presence of flammable gas.
4. Do **NOT** drop this unit onto a hard surface.
5. Do **NOT** attempt to dismantle the Obi TENS.
7. Only use specified batteries and electrodes.
8. If damaged, do not use. Return to supplier.
9. Remove the batteries when not in use.
10. Do **NOT** use while driving or operating potentially dangerous machinery or while using a microwave.
11. Keep out of the reach of children.
12. Do **NOT** place electrodes on or near the eyes, in the mouth, over the front or sides of the neck, across the head, heart, or an area of broken, inflamed, infected or numb skin. Electrodes should generally only be applied to skin with normal sensation unless under medical supervision.
13. Do **NOT** use in the presence of tuberculosis, malignant tumours, very high or very low blood pressure, high fever or acute inflammatory disease unless under medical supervision.
14. A mild shock may be caused if one electrode is accidentally removed. In the event of this happening the unit should be turned off immediately.

- b) Now insert the red and white pins (at the other end of the leadwire) into the sockets attached to all of the sticky electrodes. You will see that the lead wire is made up of 4 wires; 2 mauve and 2 white. Do not remove the electrodes from the backing at this stage (see figure). NB. Each leadwire has a red and white pin at the end. It makes NO difference whatsoever which electrode they are attached to!

- c) If you wish, attach the easy-release neck cord which enables you to hang the Obi TENS from your neck whilst in labour.



Insert a red and white pin (at the end of the lead wire) into each pair of electrodes.